



<https://en.bmc-medical.com>



G3 C20 | G3 A20

Intelligent device knows you better



CPAP / APAP


The only Chinese company with FDA approval and active sales in the United States.




Trusted Performance Tailored Comfort

Higher Treatment Efficacy with Lower Pressure Burden

- Delivers market-leading efficacy in correcting obstructive apneas and hypopneas.
- Achieves clinical targets with lower peak pressures, improving comfort and adherence.
- Precisely distinguishes extensive respiratory events to make corresponding pressure strategies, resolving unique issues for each individual.

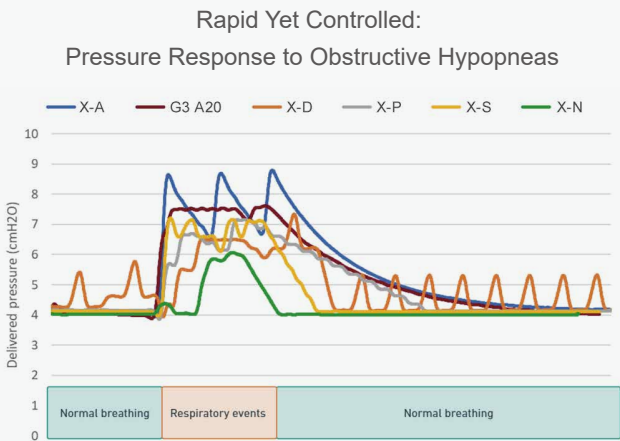


G3 stands out by delivering a prompt but measured response to obstructive hypopneas without overshooting—correcting events effectively while maintaining patient comfort.

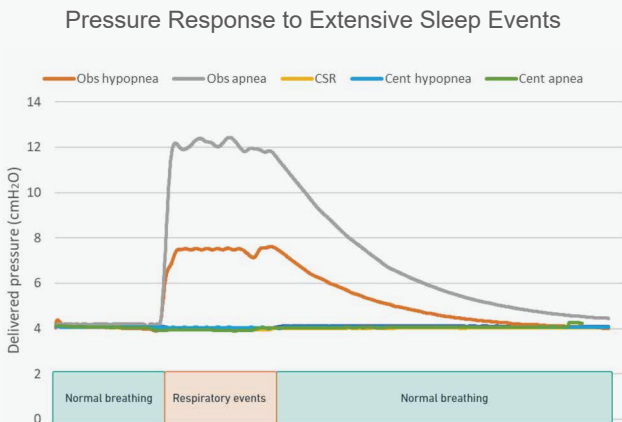

G3 accurately detects a wide range of respiratory events and delivers targeted pressure adjustments exclusively for obstructive events.

Industry-Leading Accuracy in Event Detection

- Enables targeted therapy with 100% detection of central apneas, hypopneas, and Cheyne-Stokes respiration (CSR), validated in independent third-party testing.
- Maintains therapy performance with a consistently high event correction rate, even during leaks.
- Blocks inappropriate pressure increases during central events, avoiding sleep disruptions and treatment emergent CSA.
- Accurately classifies mild/moderate/severe flow limitations with graded pressure responses, preventing progression into apneas and hypopneas.

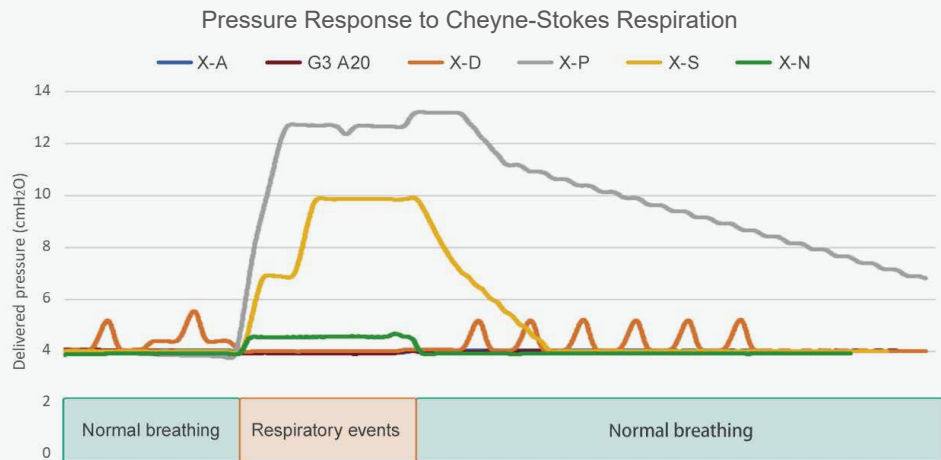



G3 correctly flags Cheyne-Stokes respiration (CSR) and prevents unnecessary pressure increases.

*Source: Data on file. Bench testing performed by KerNel Biomedical, Rouen, France (2025). Evaluation compared auto-adjusting PAP devices and assessed their responsiveness to simulated sleep events on a test bench.

Your Personalized Sleep Therapy


Hidden Doctor Assistant - Smart Adaptive Algorithm

- Dynamically adapts to an ideal starting pressure after ramp, based on data from your last 5 effective treatment days, reducing the risk of residual events emerging at the start of treatment
- Continuously optimizes minimum and maximum pressure thresholds throughout long-term therapy to accommodate weight or health fluctuations, reducing clinical visits while maintaining treatment efficacy

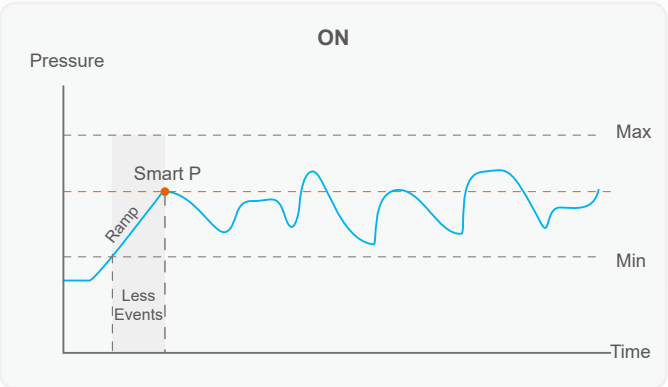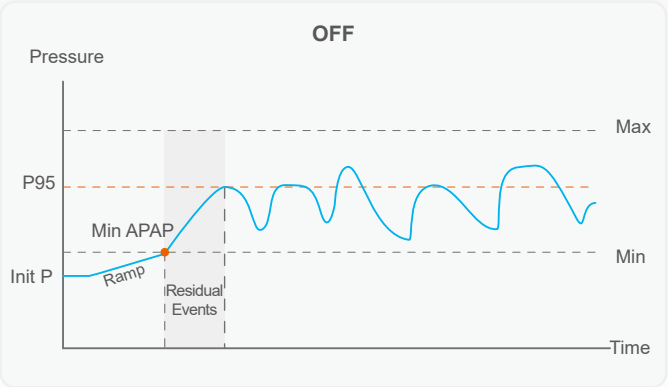


Matching different user sensitivities - 3-Level Pressure Response

- 3-Level pressure response (Fast/Standard/Soft) adapts to all OSA severity levels and pressure tolerances
- The SOFT option is ideal for pressure-sensitive patients and women—its gentle, gradual pressure transitions minimize discomfort to enhance therapy adherence.
- The FAST option is designed for patients with severe OSA and elevated BMI—delivering rapid pressure response to arrest respiratory events for higher efficacy

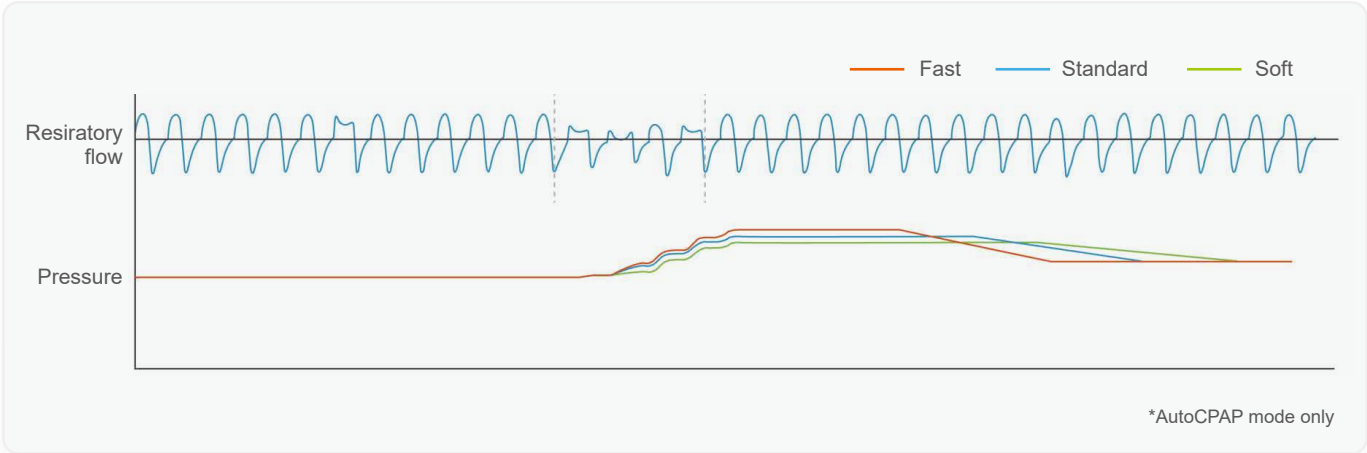


Warm to Soft Automated Care from Dusk to Dawn

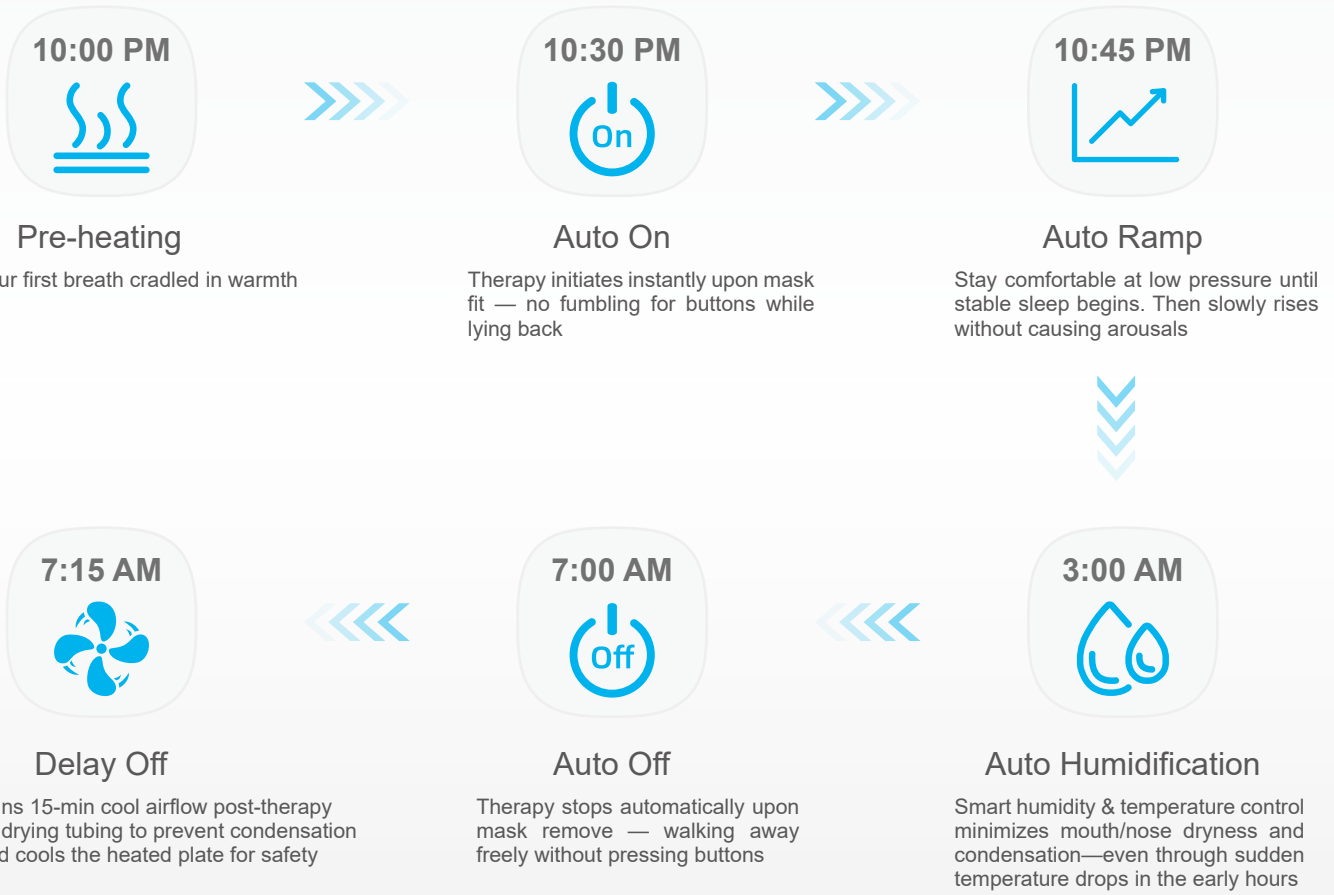

Strong Humidification All Night Moisture

360mL Oversized Chamber:

Enables 10-hour consistent humidity during therapy — no refills, no decline in comfort

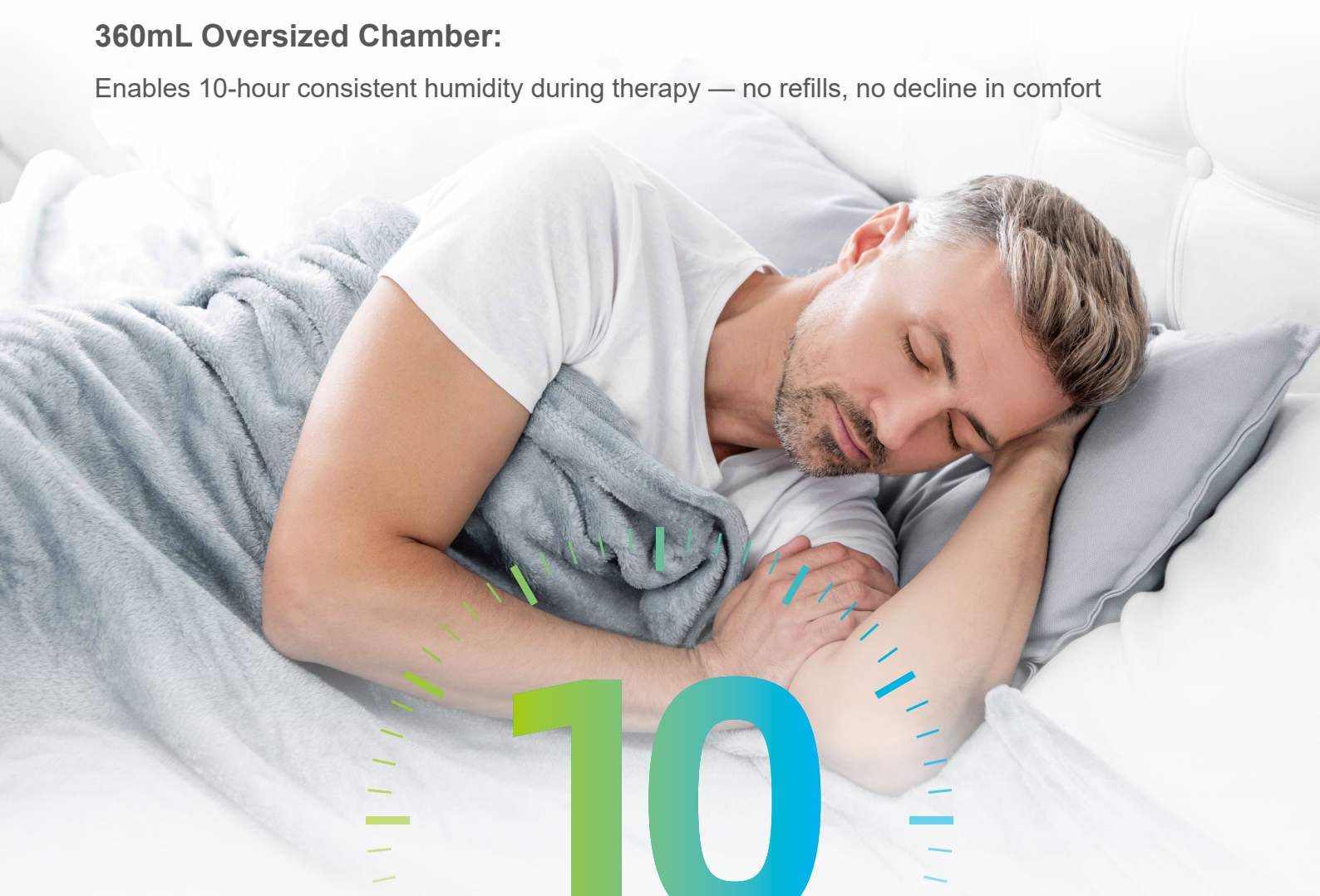

10
hours



Water Chamber with Large Capacity



Smart Control of Humidity & Temperature

Last for 10 hours

*Tested with Room temperature: 25°C, Ambient humidity: <20%, Humidity setting: Level 3

Push-Push Water Chamber Effortless Daily Refills

- **Easy remove with a simple palm press**

Ideal for users with limited finger strength



- **Compact bedside fit**

Lift-up removal saves side space vs. traditional pull-out designs



- **Extra-wide fill port**

Spill-resistant design for mess-free refilling



Data-Driven Cloud Management From Compliance to Confidence

PAP Link Digital Health solution, comprising Web, PC and APP versions, provides you with a complete set of tools to increase efficiencies by streamlining workflows, drive better patient outcomes through systematic management, and boost their therapy comfort.

For Clinicians

Precision Treatment at Your Fingertips

Remotely fine-tune settings during tricky first weeks, delivering personalized therapy to accelerate comfort and build confidence

For Distributors

Data-Driven Patient Advocacy

Cloud analytics spotlight compliance trends, enabling proactive interventions for at-risk patients;
Automated reimbursement-ready reports streamline paperwork with one-click exports;

For Patients

Therapy Insights, Anytime Anywhere

Nightly reports with smart interpretations and personalized coaching tips, empowering self-management for better sleep



PAP Link App



Download App

2015

Started from year

100+

Countries covered


1,300,000+

Registered patients






Accessories




AF-D1 Black

Air Filter




Black

Air Filter Cover




PF-D1 Black

PM 2.5 Filter




LH1 19mm

Heated tubing




SP-100

**SpO₂ kit
(Main Device)**




SOO10B-S

**Oximeter Sensor
(Finger Clip)**




SOO26N-L

**Oximeter Sensor
(Disposable)**




WT-D1

Water Chamber




CC-D1


Carrying Case



Power Adapter



Power Cable



SD Card

Specifications

	G3 C20	G3 A20
Work Mode	CPAP	CPAP, AutoCPAP
Pressure Range	4 to 20 hPa	4 to 20 hPa
Prescription Guideline		
Smart Adaptive	SmartC	SmartC, SmartA
Pressure Response	/	Fast, Standard, Soft
Comfort Functions		Device Information
Pre-heating	Support; 30 min countdown (can stop at any time)	Dimensions265 mm × 145 mm × 114 mm (L x W x H)
Ramp	Auto, 0-60 minutes	Weight1.7 kg
Reslex	Off, Level 1 - 3	Sound Pressure Level28±2 dB(A)
Auto On/Off	Support	Water capacity360±30 mL
Humidification	Auto, Level 1 - 5; Heated tubing (Optional)	Humidifier Output>12 mg H ₂ O/L
Data Communication		
Data Transfer	SD Card, iCode QR, iCode QR+, Cellular Module, WiFi Module	
Internal Storage	Summary data: up to 512 days	Detailed data: up to 12 hours
SD Card Storage	Summary data: all the time	Detailed data: up to 70 days (at 8 hours of daily use)