



## Full Face Mask (F6)

To redefine the comfort of full face therapy  
family cares



## Main features



### Under-the-nose Cushion

The under-the-nose cushion brings liberty for the nasal bridge and meanwhile minimizes the claustrophobia.



### Circular Dispersion Vent

The circular dispersion vent keep the noise as low as approx. 17dB (A) and meanwhile protect your bed partner from being disturbed.



### Modular Main Frame

The main frame is compatible with all cushion sizes (S, M, L, MW) of original packing.



### Crown-style Adjustable Headgear

The adjustable headgear is fit for most of the head sizes.




### 360° Rotatable Elbow

More freedom

## Technical Specifications

Pressure-Flow  
Curve

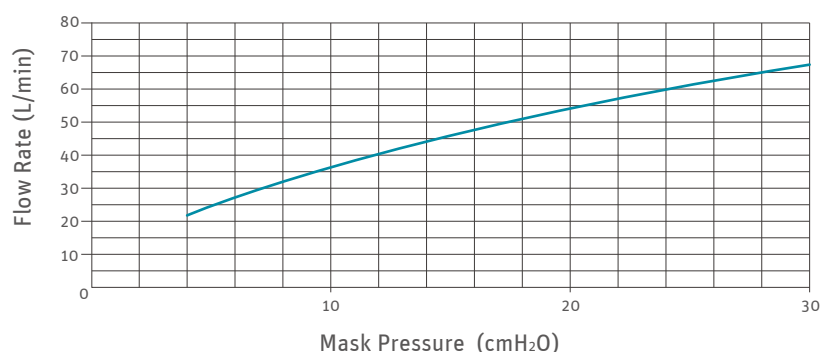

Pressure (cmH <sub>2</sub> O)	4	10	12	17	24	30
Flow Rate (L/min)	22	37	41	49	59	67

### Dead Space Information

Dead space is the empty volume of the mask up to the swivel. The dead space of the mask varies according to cushion sizes but is less than 180ml.

### Therapy Pressure

4 to 30 cmH<sub>2</sub>O

### Resistance

Drop in Pressure measured  
at 50 L/min: 0.15 cmH<sub>2</sub>O, at 100 L/min: 0.5 cmH<sub>2</sub>O

### Sound

DECLARED DUAL-NUMBER NOISE EMISSION VALUES in accordance with ISO 4871. The A-weighted sound power level of the mask is 25 dBA, with uncertainty 3 dBA. The A-weighted sound pressure level of the mask at a distance of 1 m is 17 dBA, with uncertainty 3 dBA.

## BMC Medical Co., Ltd.

Room 110, Tower A Fengyu Building, No.115 Fucheng Road, Haidian, 100036 Beijing, PEOPLE'S REPUBLIC OF CHINA

Tel: +86-10-5166 3880 / Fax: +86-10-5166 3880 ext. 810 / [www.bmc-medical.com](http://www.bmc-medical.com) / V1.1 SC-A-202307011

The product may NOT be available in your market. Please contact us for further information.