



Nasal Mask (N6)

To remold the image of nasal mask

family cares



Main features




Under-the-nose Cushion

The innovative under-the-nose cushion minimizes the contact with your nose so that no pressure needs to be exerted on the nose bridge.




Circular Dispersion Vent

The circular dispersion vent keep the noise as low as 18 dB (A) and meanwhile protect your bed partner from being disturbed.

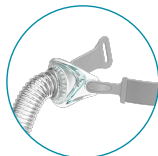

Modular Main Frame


The main frame is compatible with all cushion sizes of original packing.

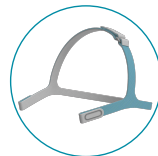

Three-point-adjustable Headgear

The adjustable headgear is fit for most of the head sizes. Wider range, more inclusiveness.

Technical specifications

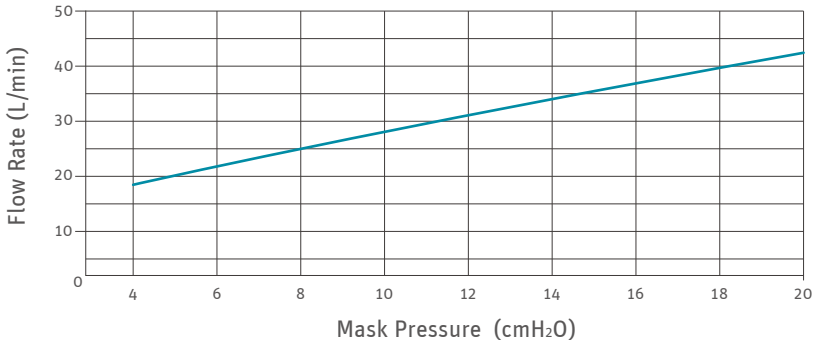Pressure-Flow Curve	<div><table><tr><th>Pressure (cmH₂O)</th><td>4</td><td>8</td><td>10</td><td>12</td><td>16</td><td>18</td><td>20</td></tr><tr><th>Flow Rate (L/min)</th><td>17</td><td>25</td><td>28</td><td>31</td><td>37</td><td>40</td><td>42</td></tr></table></div>	Pressure (cmH ₂ O)	4	8	10	12	16	18	20	Flow Rate (L/min)	17	25	28	31	37	40	42
Pressure (cmH ₂ O)	4	8	10	12	16	18	20										
Flow Rate (L/min)	17	25	28	31	37	40	42										
Therapy Pressure	4 to 20 cmH ₂ O.																
Resistance	Drop in Pressure measured (nominal) at 50 L/min: 0.4 cmH ₂ O at 100 L/min: 1.2 cmH ₂ O																
Sound	DECLARED DUAL-NUMBER NOISE EMISSION VALUES in accordance with ISO 4871. The A-weighted sound power level of the mask is 26 dBA, with uncertainty 3 dBA. The A-weighted sound pressure level of the mask at a distance of 1 m is 18 dBA, with uncertainty 3 dBA.																

BMC Medical Co., Ltd.

Room 110, Tower A Fengyu Building, No.115 Fucheng Road, Haidian, 100036 Beijing, PEOPLE'S REPUBLIC OF CHINA

Tel: +86-10-5166 3880 / Fax: +86-10-5166 3880 ext. 810 / www.bmc-medical.com / V1.1 SC-A-202306004

Please contact us for more details, such as the available models and accessible market.